

AUTO hMeDIP KIT MANUAL

Hydroxymethylated DNA Immunoprecipitation

Auto hMeDIP kit x16 (*monoclonal rat antibody*)

Cat. No. C02010033 (*Old: AF-Auto02-0016*)

Auto hMeDIP kit x16 (*monoclonal mouse antibody*)

Cat. No. C02010034

Auto hMeDIP kit x16 (*polyclonal rabbit antibody*)

Cat. No. C02010035

FEATURES

- Get your hydroxymethylation profile in 24h
- Includes hydroxymethylated, methylated and unmethylated DNA for QC
- Improved reproducibility using the IP-Star®
- Highly specific (monoclonal antibody for 5hmc)
- Optimized DNA isolation buffer

Technical Assistance & Ordering Information

Diagenode s.a. BELGIUM | EUROPE

LIEGE SCIENCE PARK
Rue Bois Saint-Jean, 3
4102 Seraing - Belgium
Tel: +32 4 364 20 50
Fax: +32 4 364 20 51
techsupport@diagenode.com
orders@diagenode.com

Diagenode Inc. USA | NORTH AMERICA

400 Morris Avenue, Suite #101
Denville, NJ 07834 - USA
Tel: +1 862 209-4680
Fax: +1 862 209-4681
techsupport.na@diagenode.com
orders.na@diagenode.com

For a complete listing of Diagenode's international distributors, visit:
<http://www.diagenode.com/en/company/distributors.php>
For the rest of the world, please contact Diagenode s.a.

Technical Assistance

At DIAGENODE we pride ourselves on the quality and availability of our technical support. Our Technical Services Departments are staffed by experienced scientists with extensive practical and theoretical expertise in molecular biology and the use of DIAGENODE products. If you have any questions, or experience any difficulties regarding the SX-8G IP-Star or DIAGENODE products in general, do not hesitate to contact us.

DIAGENODE customers are a major source of information regarding advanced or specialized uses of our products. This information is helpful to other scientists as well as to the researchers at DIAGENODE. We therefore encourage you to contact us if you have any suggestions about product performance or new applications and techniques.

For technical assistance and more information call the DIAGENODE Technical Service Department or contact your local distributor.

Contents

- Introduction..... 4**
- SX-8G IP-Star Automated System for ChIP, MeDIP &MBD..... 5**
- Kit Method Overview 8**
- Kit Materials 9**
 - Kit Content 9
- How to perform Automated IPure in the SX-8G IP-Star® Compact..... 11**
 - Running a protocol..... 12
- How to perform Automated hMeDIP in the SX-8G IP-Star® 18**
 - Loading and running protocol..... 19
 - Shutting down the SX-8G IP-Star® 22
- Quantitative PCR & Data Analysis 23**
- Results..... 24**
- Troubleshooting Guide..... 24**
- Ordering Information..... Back Cover**

Introduction

The Novel Diagenode Auto hMeDIP kit is designed to immunoprecipitate hydroxymethylated DNA (hMethyl DNA IP). This kit is the first and only example of MeDIP kit specifically designed and fully validated for affinity-capture and detection of hydroxymethylated regions using the Diagenode highly specific antibody against 5-hmc.

One of the fastest growing fields in biology and cancer research is epigenetics. While the underlying genetic code defines which proteins and gene products are synthesized, it is epigenetic control that defines when and where they are expressed. Epigenetic control is generally mediated by methylation of cytosine to 5-methylC (5-mC) in CpG islands and post-translational modification of histones. Methylation of CpGs near promoters is associated with gene silencing, as is deacetylation of histones.

There is substantial interest and speculation in the role of a recently discovered second type of DNA methylation, 5-hydroxymethylcytosine (5-hmC), although its precise function has not yet been elucidated. This new cytosine base modification results from the enzymatic conversion of 5-Methylcytosine into 5-Hydroxymethylcytosine by the TET family of oxygenases. Preliminary results indicate that 5-hmC may have important roles distinct from 5-mC. Although its precise role has still to be shown, early evidence suggests a few putative mechanisms that could have big implications in epigenetics: 5-hydroxymethylcytosine may well represent a new pathway to demethylate DNA involving a repair mechanism converting hmC to C and, as such open up entirely new perspectives in epigenetic studies.

Due to the structural similarity between 5-mC and 5-hmC, these bases are experimentally almost indistinguishable. Recent articles demonstrated that the most common approaches (eg. enzymatic approaches, bisulfite sequencing) do not account for 5-hmC. The development of the affinity-based technologies appears to be the most powerful and so far the only way to differentially and specifically enrich 5mC and 5hmC sequences.

Recently, Diagenode launched new highly-specific monoclonal antibodies and kits for the differential study of the functions of 5-hmC and 5-mC. In the Auto hMeDIP kit, our antibody directed against 5-hydroxymethylcytosine is provided as well as hmeDNA, meDNA and unDNA internal IP PCR controls. The IP has been optimized to specifically select and precipitate hydroxymethylated DNA fragments by the use of our antibody, buffers and protocol. The IP efficiency can indeed be doublechecked with the use of our internal controls.

The Auto hMeDIP kits allow you to perform DNA hydroxymethylation analysis of your sample together with optimized internal IP control. Performing hydroxymethylation profiling with the Auto hMeDIP kit is FAST, RELIABLE and HIGHLY SPECIFIC.

In the Auto hMeDIP kit, the protocol has been improved to allow researchers to work in small tubes. The kit ensures the use of low amount of reagents per reaction. The number of steps is reduced and handling is easier with our Magnetic hMethyl DNA IP procedure. The hMeDIP protocol is flexible, as the IP'd DNA can be isolated and/or purified in different ways based on the downstream application: an extra-fast and simplified protocol is included in the Auto hMeDIP kit (for qPCR analysis); but more traditional methods are also proposed (see separate module). Moreover, our Magnetic Rack together with our new hMeDIP kit protocol ensures the best IP conditions by working at a constantly cooler temperature. The Diagenode Magnetic Rack has been designed to be used in IP experiments, keeping samples cool longer and allowing the use of small tubes to reduce the reaction volumes and waste of reagents.

The hMeDIP kit protocol has been validated with our Bioruptor®. Nevertheless, DNA can be sheared with any in house protocol and sonication apparatus as long as efficiency is checked before use.

The Auto hMeDIP kit contains all reagents you need for you hMeDIP Assay but it can be purchased with two additional modules for: 1/ Preparation of larger quantities of genomic DNA and 2/ "Traditional" purification of the IP'd DNA for subsequent next generation sequencing analysis.

1. The XL GenDNA Module is optimized for the preparation of large quantities of DNA ready-to-use in MeDIP. An optimized protocol for DNA shearing is provided as well.
2. The DNA purification includes all the reagents and buffers needed: eluting DNA from the washed beads after IP (using buffers D, E, F) and to proceed to phenol/chloroform extractions and ethanol precipitation (using DNAIP TE, DNA-IP co-precipitant and DNA-IP precipitant). Note that purification columns can also be used after elution.

SX-8G IP-Star® and SX-8G IP-Star® Compact Systems for automation of epigenetic applications

Diagenode has developed two automated platforms (SX-8G IP-Star® and SX-8G IP-Star® Compact) designed to increase your lab's productivity, efficiency and experimental reproducibility. The two automated platforms are capable of processing up to 16 samples per cycle. The automated systems process sheared chromatin (or DNA) to deliver purified DNA ready for qPCR, amplification, microarray and sequencing analysis. Both, the SX-8G IP-Star® and SX-8G IP-Star® Compact have an easy-to-use open software that provides you with flexibility. This allows you to create your personal protocol according to your specific needs.

Major benefits of Diagenode Automated Platforms



SX-8G IP-Star® Compact



SX-8G IP-Star®



- High resolution ChIP-seq and MeDIP-seq profiles
- Automated library preparation for Next Generation sequencing
- Reduces hands on time to just 30 minutes
- Reduces variability between operators and labs
- Ideal for low sample starting amounts
- Compatible with Diagenode Kits (Auto ChIP kit, Auto Histone ChIP-seq kit, Auto Transcription ChIP-seq kit, Auto True MicroChIP kit, Auto MeDIP kit, Auto MethylCap kit, Auto hMeDIP, Auto IPure kit)
- Reduces cross-contamination

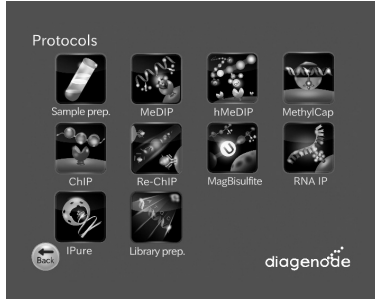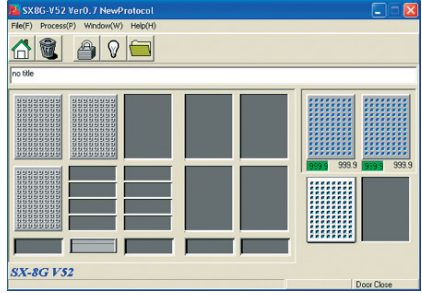	SX-8G IP-Star® Compact	SX-8G IP-Star®
Applications	ChIP-seq, MeDIP-seq, MethylCap-seq, hMeDIP, IPure, Sample preparation, Re-ChIP, MagBisulfite, RNA-IP, Library preparation for NGS platforms.	ChIP-seq, MeDIP-seq, MethylCap-seq, hMeDIP, IPure, Sample preparation, Re-ChIP, MagBisulfite, RNA-IP.
Software		
User interface	Intuitive touch screen panel	PC Software
User friendly	Software training not required	Software training before use
Dispensing	Automated dispensing of assay reagents	Manual dispensing of assay reagents
Protocol optimization (flexible parameters)	Antibody coating (temperature, time, mixing speed) Immunoprecipitation (temperature, time, mixing speed) Washes (temperature, time, mixing speed)	Antibody coating (temperature, time) Immunoprecipitation (temperature, time)
New protocol development	Achievable by Diagenode product specialist	Achievable by customer after training
Characteristics	750W x 740 D x 610 H 100 kg 8 Nozzles X-Y-Z axis 4 – 95°C	1070W x 650 D x 780 H 130 kg 8 Nozzles X-Y-Z axis 4-95°C

Improved reproducibility

Our SX-8G IP-Star will increase the immunoprecipitation reproducibility between IPs performed by the same as well as by different operators (see figure 1 and 2 below). Reagents (Antibodies, buffers,...) and sheared chromatin were identical for "ManChIP" and "AutoChIP". The SX-8G IP-Star Automated system removes variation that can be created by manual handling and allows you to optimize and standardize your assay within a lab. The SX-8G IP-Star is designed to improve the accuracy and the reproducibility of any immunoprecipitation experiment.

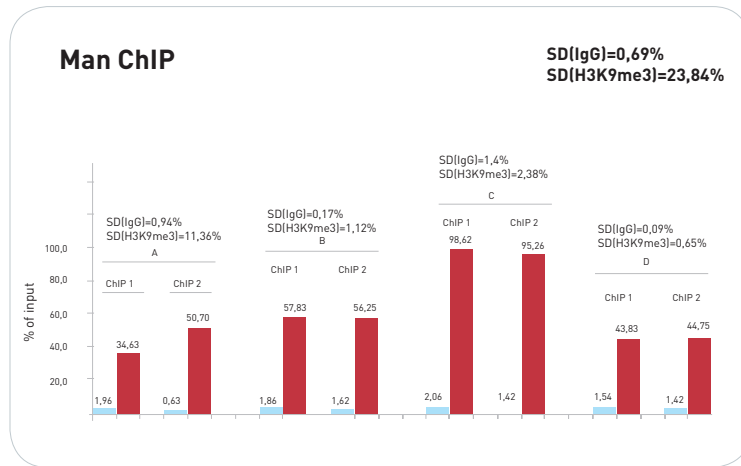


Figure 1: Manual ChIP. Four different operators have each performed two ChIP experiments using H3K9me3 antibody on the genomic region SAT2 (positive locus). 10,000 HeLa cells have been used per IP. Reagents and sheared chromatin were identical per assay. The standard deviations between the ChIPs performed by the same operator and between the four different operators are displayed.

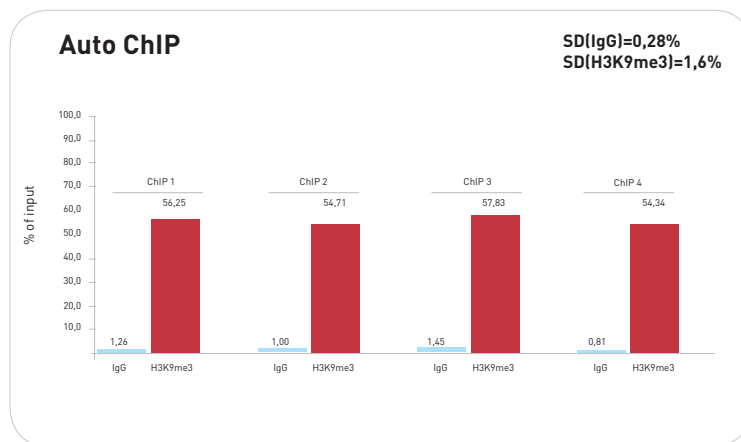


Figure 2: Automated ChIP. Four ChIP experiments using H3K9me3 antibody on the genomic region SAT2 (positive locus) have been performed by the SX-8G IP-Star. 10,000 HeLa cells have been used per IP. Reagents and sheared chromatin were identical per assay. The standard deviations between the four ChIPs performed by the SX-8G IP-Star are displayed.

Kit Method Overview

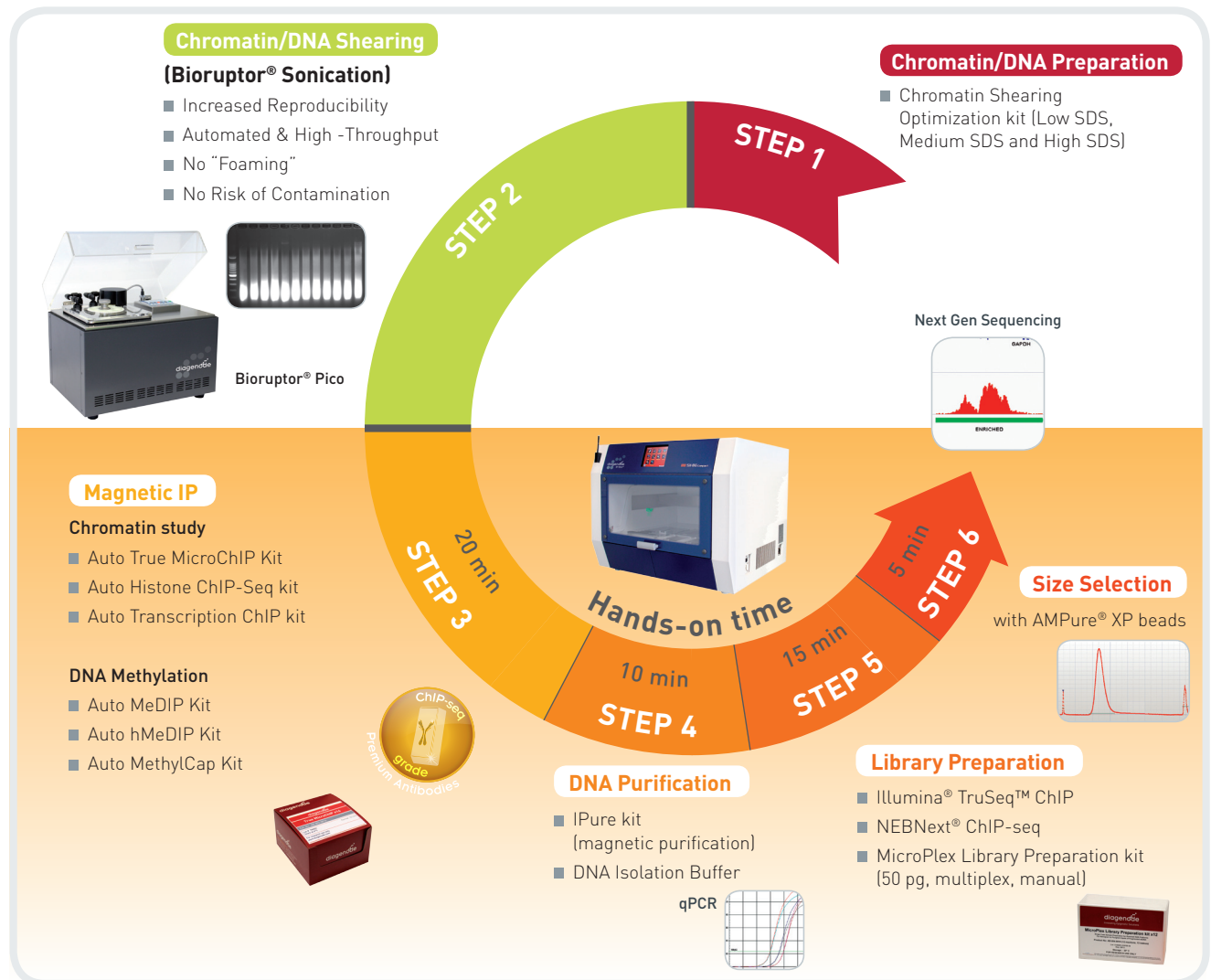


Figure 3. Diagenode provides a full suite of automated solutions for ChIP experiments.

For Step 1, we offer products to isolate nuclei and chromatin. Step 2 describes reproducible sample shearing with the Bioruptor® product line. In Step 3 and Step 4, the Diagenode IP-Star Compact provides error-free, walk-away automation for all your immunoprecipitation and antibody capture needs.

Kit Materials

Kit Content

Auto hMeDIP kit x16 (monoclonal rat antibody) Cat. No. C02010033 (Old: AF-Auto02-0016)				
Component	Cat. No.	Comments	Quantity	Storage
Water	EF-113-0008	-	24 ml	4°C
hMeDIP buffer H1	EF-114-0001	Ion chelator mix included (10x).	3 ml	4°C
hmeDNA control	EF-115-0040	PCR product (0,02 ng/μl).	40 μl	-20°C
meDNA control	EF-116-0040	PCR product (0,02 ng/μl).	40 μl	-20°C
unDNA control	EF-117-0040	PCR product (0,02 ng/μl).	40 μl	-20°C
5-hmC monoclonal antibody (rat)	EF-118-0032	1,6 μg/μl	32 μl (50 μg)	-20°C
Rat IgG	EF-119-0050	1 μg/μl	50 μl	-20°C/4°C
Protein G-coated magnetic beads	EF-120-0220	The beads are supplied for 16 IPs; detergent and 0.02%, sodium azide included.	220 μl	4°C Do not freeze
hMeDIP buffer H2	EF-121-0004	BSA and Ion chelator mix included.	12 ml	4°C
hMeDIP buffer H3	EF-122-0002	Ion chelator mix included.	4 ml	4°C
DNA isolation Buffer (DIB)	EF-123-0004	-	5 ml	4°C
Proteinase K	EF-124-0040	100 x stock solution.	40 μl	- 20°C
hmeDNA primer pairs	EF-125-0050	5 μM each (Rv & Fw).	50 μl	- 20°C
meDNA primer pairs	EF-126-0050	5 μM each (Rv & Fw).	50 μl	-20°C
unDNA primer pairs	EF-127-0050	5 μM each (Rv & Fw).	50 μl	-20°C
Mouse Sfi1 primer pairs	EF-128-0050	5 μM each (Rv & Fw).	50 μl	-20°C

Auto hMeDIP kit x16 (monoclonal mouse antibody) Cat. No. C02010034				
Component	Cat. No.	Comments	Quantity	Storage
Water	EF-113-0008	-	24 ml	4°C
hMeDIP buffer H1	EF-114-0001	Ion chelator mix included (10x).	3 ml	4°C
hmeDNA control	EF-115-0040	PCR product (0,02 ng/μl).	40 μl	-20°C
meDNA control	EF-116-0040	PCR product (0,02 ng/μl).	40 μl	-20°C
unDNA control	EF-117-0040	PCR product (0,02 ng/μl).	40 μl	-20°C
5-hmC monoclonal antibody (mouse)	EF-131-0050	1 μg/μl	50 μl (50 μg)	-20°C
Mouse IgG	EF-119-0050	1 μg/μl	50 μl	-20°C/4°C
Anti-mouse IgG-coated magnetic beads	EF-133-0220	The beads are supplied for 16 IPs; detergent and 0.02%, sodium azide included.	220 μl	4°C Do not freeze
hMeDIP buffer H2	EF-121-0004	BSA and Ion chelator mix included.	12 ml	4°C
hMeDIP buffer H3	EF-122-0002	Ion chelator mix included.	4 ml	4°C
DNA isolation Buffer (DIB)	EF-123-0004	-	5 ml	4°C
Proteinase K	EF-124-0040	100 x stock solution.	40 μl	- 20°C
hmeDNA primer pairs	EF-125-0050	5 μM each (Rv & Fw).	50 μl	- 20°C
meDNA primer pairs	EF-126-0050	5 μM each (Rv & Fw).	50 μl	-20°C
unDNA primer pairs	EF-127-0050	5 μM each (Rv & Fw).	50 μl	-20°C
Mouse Sfi1 primer pairs	EF-128-0050	5 μM each (Rv & Fw).	50 μl	-20°C

Auto hMeDIP kit x16 (polyclonal rabbit antibody) Cat. No. C02010035				
Component	Cat. No.	Comments	Quantity	Storage
Water	EF-113-0008	-	24 ml	4°C
hMeDIP buffer H1	EF-114-0001	Ion chelator mix included (10x).	3 ml	4°C
hmeDNA control	EF-115-0040	PCR product (0,02 ng/μl).	40 μl	-20°C
meDNA control	EF-116-0040	PCR product (0,02 ng/μl).	40 μl	-20°C
unDNA control	EF-117-0040	PCR product (0,02 ng/μl).	40 μl	-20°C
5-hmC polyclonal antibody (rabbit)	EF-134-0050	crude serum	50 μl	-20°C
Rabbit IgG	EF-119-0050	1 μg/μl	50 μl	-20°C/4°C
Protein A-coated magnetic beads	EF-136-0220	The beads are supplied for 16 IPs; detergent and 0.02%, sodium azide included.	220 μl	4°C Do not freeze
hMeDIP buffer H2	EF-121-0004	BSA and Ion chelator mix included.	12 ml	4°C
hMeDIP buffer H3	EF-122-0002	Ion chelator mix included.	4 ml	4°C
DNA isolation Buffer (DIB)	EF-123-0004	-	5 ml	4°C
Proteinase K	EF-124-0040	100 x stock solution.	40 μl	- 20°C
hmeDNA primer pairs	EF-125-0050	5 μM each (Rv & Fw).	50 μl	- 20°C
meDNA primer pairs	EF-126-0050	5 μM each (Rv & Fw).	50 μl	-20°C
unDNA primer pairs	EF-127-0050	5 μM each (Rv & Fw).	50 μl	-20°C
Mouse Sfi1 primer pairs	EF-128-0050	5 μM each (Rv & Fw).	50 μl	-20°C

Components available separately			
Description	Reference	Quantity	Storage
96 well microplates	WA-003-0010	10	RT
Tips (bulk)	WC-001-1000	1000	RT
200 μl tube strips (12 tubes/strip) + cap strips	WA-001-0080	80	RT
200 μl tube strips (8 tubes/strip) + cap strips for SX-8G IP-Star® Compact	WA-002-0120	120	RT
Tips (box)	WA-002-960	10 x 96	RT

Kits and Modules available separately		
Description	Reference	Quantity
Chromatin shearing optimization kit - Low SDS	AA-001-0100	1 kit
Chromatin shearing optimization kit - Medium SDS	AA-002-0100	1 kit
Chromatin shearing optimization kit - High SDS	AA-003-0100	1 kit

IPure		
Description	Reference	Quantity
IPure	AL-100-0100	100 rxns
Auto IPure	AL-Auto01-0100	100 rxns
96 well microplates	WA-003-0010	10 pc

How to perform Automated ChIP in the SX-8G IP-Star® Compact



SX-8G IP-STAR® COMPACT

How to perform Automated hMeDIP in the SX-8G IP-Star® Compact

Prior to hydroxyMethylated DNA immunoprecipitation (hydroxyMethyl DNA IP), DNA samples are first prepared and sheared with Diagenode GenDNA module (Cat. No. C03030020 (mc-magme-003)).

Prepare reagents

1. Prepare the IP incubation mix w/o antibodies and w/o magnetic beads for all your hMeDIP reactions (Table 1).


Table 1: IP incubation mix with no antibodies and no beads

Reagent	Volume per IP + INPUT	Volume per additional IP
Water	91.5 µl	76.25 µl
hMeDIP buffer H1	12 µl	10 µl
hmeDNA control	1.5 µl	1.25 µl
meDNA control	1.5 µl	1.25 µl
unDNA control	1.5 µl	1.25 µl
DNA sample (0.1 µg/µl)	12 µl	10 µl
TOTAL VOLUME	120.00 µl	100.00 µl

NOTE: hMeDIP Buffer H1 contains detergent; if the appearance is cloudy and crystallized please warm gently prior to use.

- Incubate at 95°C for 10 minutes.
 - Quickly chill sample on ice (it is best to use ice-water).
 - Perform a pulse spin to consolidate your sample.
 - First, take out 10 µl per INPUT (that is 10% input) and transfer to a new labelled tube.
 - Keep the input samples at 4°C. The input sample is to be used as a control of starting material and it is therefore not to be used in IP.
 - 100 µl of the IP incubation mix will be used per IP
2. Prepare hMeDIP H1 (1:10) by diluting hMeDIP H1 buffer in water.
300 µl hMeDIP H1 (1:10) buffer are needed per IP reaction

Running a protocol

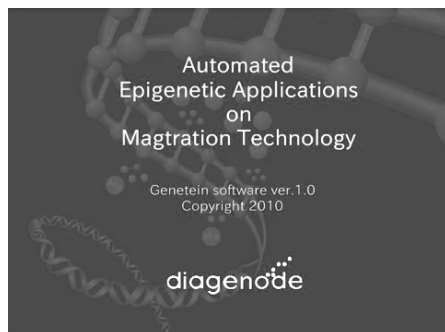

Diagenode Splash Screen – A0

After the software start-up screen disappears, the Diagenode splash screen is displayed for several seconds, and then disappears.




Start Screen – Top menu

After the Digenode splash screen disappears, the start screen is displayed. This is the first active window; it allows the user to enter into three different parts of the software.

USER ACTIONS:


Buttons:

- Protocols
- Maintenance (for technical services)
- Information (for Diagenode contact details)

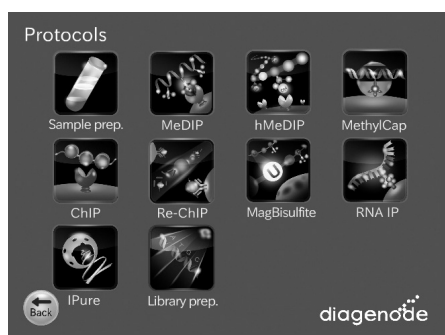

Protocols screen

All available protocols are displayed on this screen.

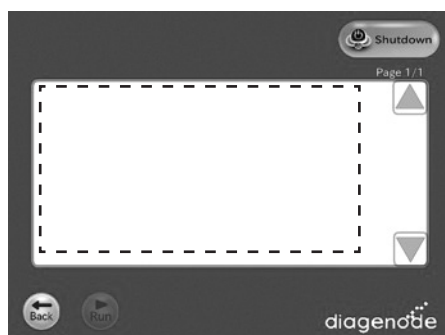


Screen – [Categories Name] Protocol List

After the user presses the "[Categories Name]" button, the "[Categories Name]" appears. When selected the protocol on the protocol list, the "Run" button shall turn executable.

Buttons:

- The user presses the "Back" button. The user returns to the "Protocols" screen.
- The user presses the "Shutdown" button. The screen shall be changed to "Power Off".
- The user presses the "Run" button. The screen shall be changed to "Sample number".
- ▲ Page up the list box.
- ▼ Page down the list box

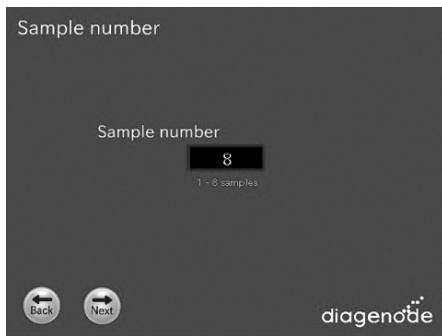


Screen – Sample number


After the user presses the “Run” button, the “Sample number” appears.

Buttons:

- The user presses the “Sample number” Text box. The screen will be changed to keyboard.
- The user presses the “Back” button. The user returns to the “Protocol List” screen.
- The user presses the “Next” button. The screen shall be changed to “Configuration” or “Layout information”.

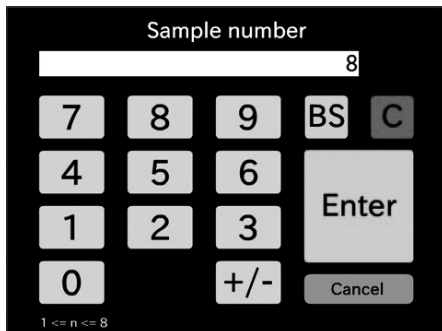

Keyboard

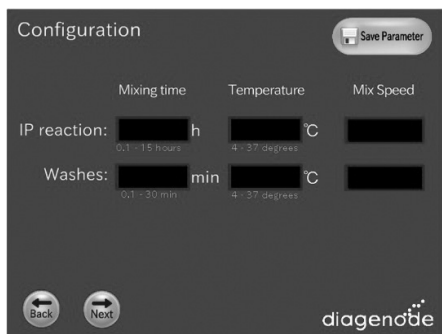


Screen – Configuration


After the user presses the next button from the “Sample number” screen, the “Configuration” screen appears.

Buttons:


- The user presses the “Back” button. The user returns to the “Protocol List” screen.
- The user presses the “Next” button. The screen shall be changed to “Layout information”.
- The user presses the “Save Parameter” button. The screen will be changed to “Save Parameter - Confirmation”.
 - OK – Current parameters shown in the Display View will be stored to the [Protocol].ptd. And, returns the user to the display of the “Configuration” screen.
 - No – Returns the user to the display of the “Configuration” screen.
- The user presses the Text box. The screen will be changed to Keyboard or Speed list menu.

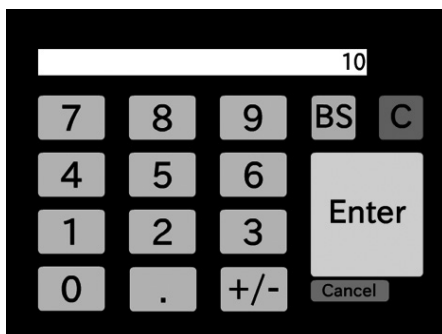


Keyboard

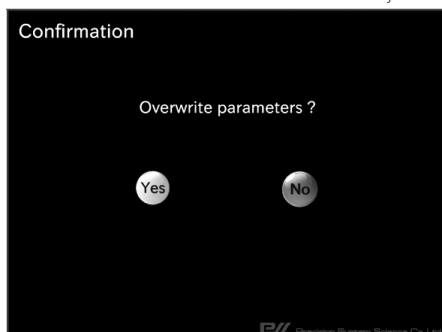


Block-PCR Tube

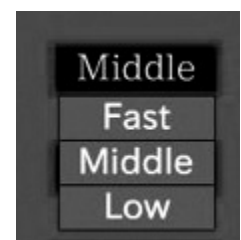


Speed list menu

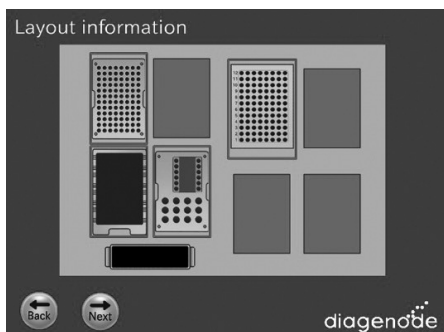


Layout information

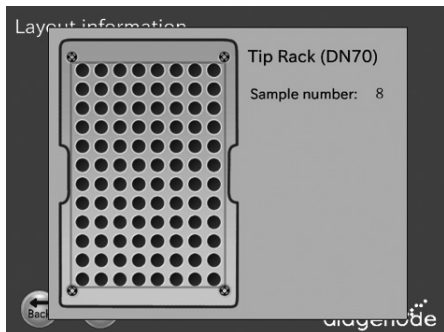


Block-Tip

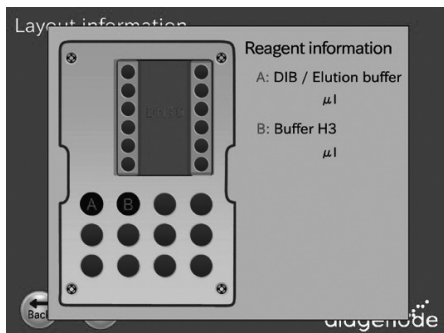

Block-Regent Tip Rack

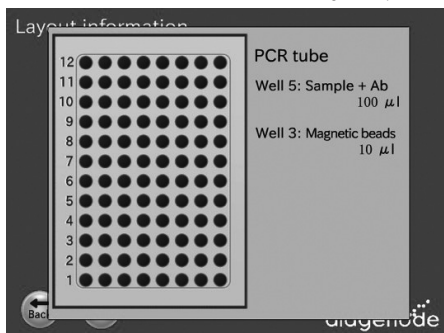

Block-PCR Tube

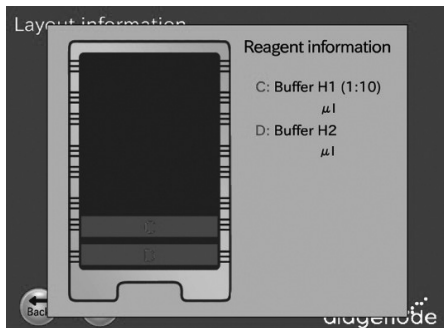

Screen – Layout Information

After the user presses the “next” button from “Sample number” screen or “Configuration” screen, the “Layout Information” screen appears.


Buttons:

- The user presses the “Back” button. The user returns to the previous screen.
- The user presses the “Next” button. The screen changes to “Set confirmation”.
- When the user presses a block, that block is magnified on the work surface layout background. The magnified view provides a better display of the correct method setup for that block on the work surface.
- Based on the selected protocol, the user follows the indications provided in the screens to set up correctly the different reagents and samples based on the selected ChIP protocol.

Screen – Layout Information

If the kit contains:

- A rat monoclonal 5-hmC antibody, add 1.6 µl of the 5-hmC antibody or 2.5 µl of Rat IgG per tube.
- A mouse monoclonal 5-hmC antibody, add 2.5 µl of the 5-hmC antibody or 2.5 µl of Rat IgG per tube.
- A rabbit polyclonal 5-hmC antibody, add 2.5 µl of the 5-hmC antibody or 2.5 µl of Rabbit IgG per tube.

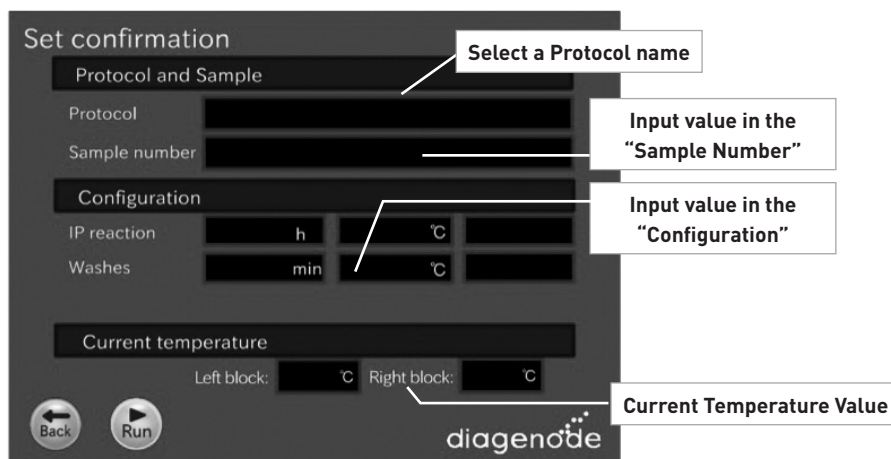


Screen – Set confirmation

After the user presses the “next” button in the “Layout information” screen, the “Set confirmation” screen appears.

At this point, user is expected to be ready to press RUN.

Buttons:

- The user presses the “Back” button. The user returns to the layout information screen.
- The user presses the “Run” button. This is the expected action when user gets to this display after reviewing blocks. Runs the protocol.

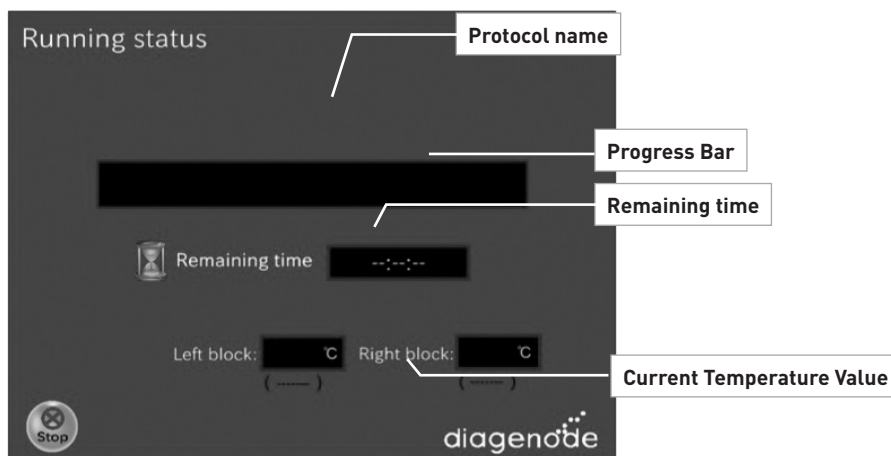


Screen – Running

After the user presses the “Run” button in the “Set confirmation” screen, the “Running” screen appears.

Buttons:

- The user presses the “Stop” button. The screen changes to “Stop Dialog”.


Status screen is preferred as a progress bar that moves across the screen as the step progresses

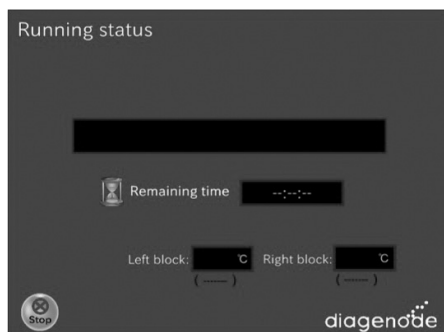

Screen – Running status

This screen gives informations about the current running step of the protocol.

The user can check through this screen the passed and remaining time of the experiment.


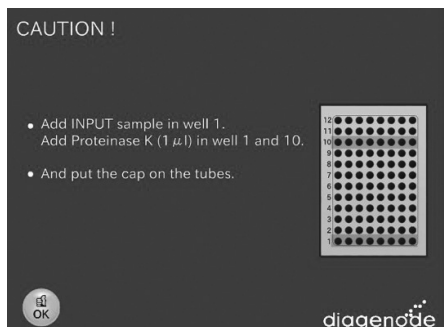


Screen – CAUTION

INPUT is defined as

INPUT= 10 μ l Incubation Mix + 90 μ l DIB buffer

IMPORTANT: Please note that the enriched methylated DNA in DIB buffer is single strand DNA that can be directly analyzed by qPCR. For downstream applications such as sequencing or arrays, the enriched methylated DNA needs to be purified by phenol/chloroform extraction and converted to double stranded DNA.

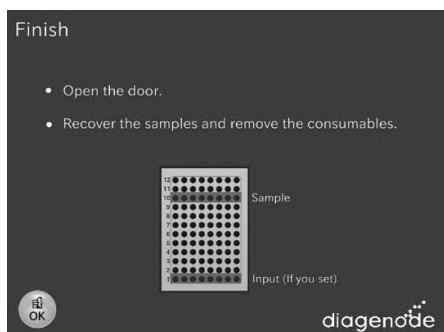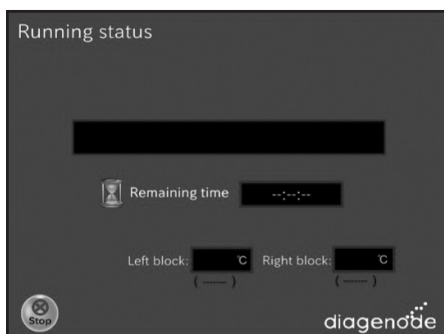

Screen – Finish/End

When the protocol is complete, a window appears telling user the run is over. The screen behind this window should be the Startup screen. When OK is pressed, then the Startup screen appears and the user can immediately begin to remove their sample and prepare the next run.



At this point, user is expected to be ready to press RUN.

Buttons:

- The user presses the "OK" button. Then screen shall be changed to "[Categories Name] Protocol List".

**Screen – Caution !**

When the protocol finishes the user can return to the protocol list (screen **A.**) or warm the peltier block (screen **B.**) to eliminate possible condensation in the block.

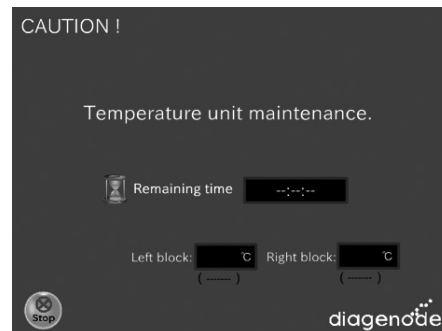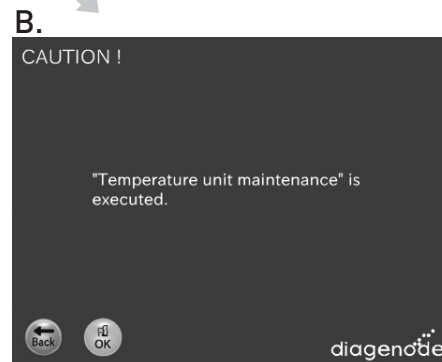

How to perform Automated ChIP in the SX-8G IP-Star®



SX-8G IP-STAR®

How to perform Automated hMeDIP in the SX-8G IP-Star®

Prior to hydroxyMethylated DNA immunoprecipitation (hydroxyMethyl DNA IP), DNA samples are first prepared and sheared with Diagenode GenDNA module (Cat# mc-magme-03).

A) Prepare reagents

1. Prepare the IP incubation mix w/o antibodies and w/o magnetic beads for all your hMeDIP reactions (Table 1).

Table 1: IP incubation mix with no antibodies and no beads


Reagent	Volume per IP + INPUT	Volume per additional IP
Water	91.5 µl	76.25 µl
hMeDIP buffer H1	12 µl	10 µl
hmeDNA control	1.5 µl	1.25 µl
meDNA control	1.5 µl	1.25 µl
unDNA control	1.5 µl	1.25 µl
DNA sample (0.1 µg/µl)	12 µl	10 µl
TOTAL VOLUME	120.00 µl	100.00 µl

NOTE: hMeDIP Buffer H1 contains detergent; if the appearance is cloudy and crystallized please warm gently prior to use.

- Incubate at 95°C for 10 minutes.
 - Quickly chill sample on ice (it is best to use ice-water).
 - Perform a pulse spin to consolidate your sample.
 - First, take out 10 µl per INPUT (that is 10% input) and transfer to a new labelled tube.
 - Keep the input samples at 4°C. The input sample is to be used as a control of starting material and it is therefore not to be used in IP.
 - 100 µl of the IP incubation mix will be used per IP
2. Prepare hMeDIP H1 (1:10) by diluting hMeDIP H1 buffer in water.
300 µl hMeDIP H1 (1:10) buffer are needed per IP reaction

B) Dispense prepared reagents into corresponding tubes (see picture below)

Loading reagents: make sure that all reagents are in the bottom of the tubes (especially magnetic beads) before starting the protocol.

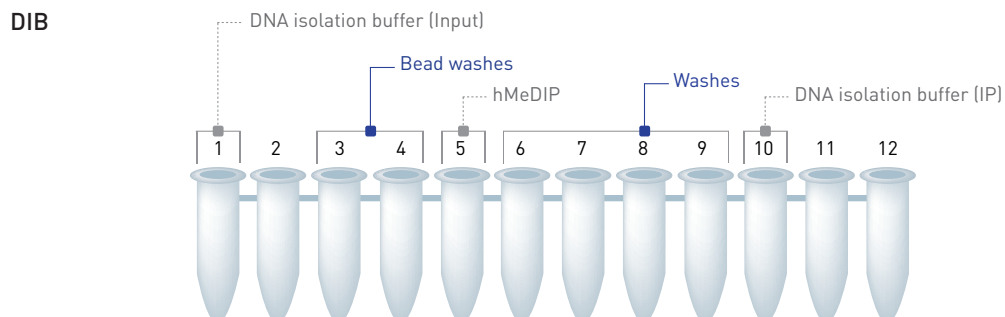

Tube #	Description	DIB
		Volume
1	DNA isolation buffer	90 µl
2	Empty	/
3	Magnetic beads	10 µl
4	hMeDIP H2	190 µl
5	IP Incubation Mix + Antibody *	100 µl
6	hMeDIP H1 (1:10)	100 µl
7	hMeDIP H1 (1:10)	100 µl
8	hMeDIP H1 (1:10)	100 µl
9	hMeDIP H3	100 µl
10	DNA isolation Buffer	100 µl
11	empty	/
12	empty	/

* If the kit contains:

- A **rat** monoclonal 5-hmC antibody, add 1.6 µl of the 5-hmC antibody or 2.5 µl of Rat IgG per tube.
- A **mouse** monoclonal 5-hmC antibody, add 2.5 µl of the 5-hmC antibody or 2.5 µl of Rat IgG per tube.
- A **rabbit** polyclonal 5-hmC antibody, add 2.5 µl of the 5-hmC antibody or 2.5 µl of Rabbit IgG per tube.

hMeDIP protocols provided for the SX-8G IP-Star

	Volumes	DIB
8 IP's	100 µl	✓
16 IP's	100 µl	✓

Loading and running protocol

Be sure that the computer connected to the robot never switches to the standby modus. (standby modus has to be inactivated). Standby of the computer will lead to the abort of the protocol.

Protocol Name	hMeDIP_14h_DIB
Reagent Preparation*	1h
Magnetic Bead Washes	2h
Immunoselection	10h
Washes	20min
Add reagents	15min
DNA isolation**	30min
Total Time	14h05 min


* Input required is sheared DNA ready-to-MeDIP

** Performed by using DIB (DNA Isolation buffer)

Note: Hands-on-work time is reduced to 1h45 min.

1. Switch on the SX-8G IP Star. The power switch is on the right side of the instrument.
2. Switch on the computer.
3. Start SX-8G V52 software through SX-8G V52 the following icon 
4. Place the prepared tube strip on the right cooling / heating block of the workstation


5. Close the workstation door and lock it using the following icon

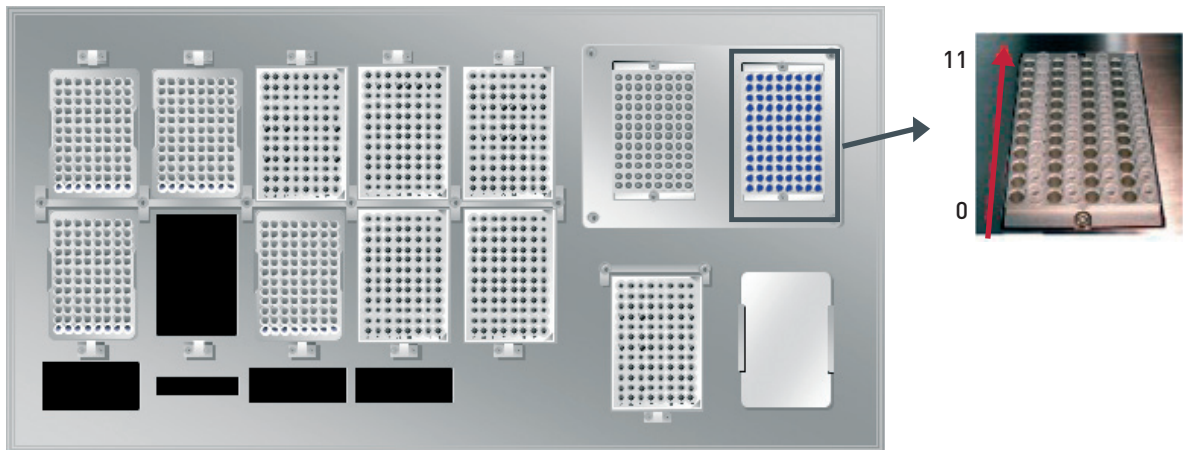

6. Press the following icon 

Select the protocol of interest. Press start. 

Before starting the protocol a start confirmation window will appear. Press OK and the protocol will run.

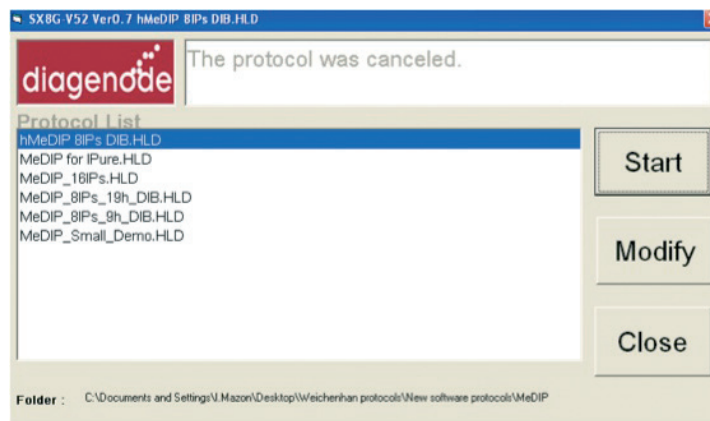

IMPORTANT NOTE:

If the ChIP protocols do not appear in the screen,


1. If the hMeDIP protocols do not appear in the screen, Open the SX-8V52 directory and open Easy start ini file. Write the directory location of the protocols
2. The Easy start ini file should contain the following information:

[EASYSTARTSCREEN]


HoldFilePath=C:\Documents and Settings\Desktop\New software protocols\ MeDIP

In red is indicated the the directory location of the MeDIP protocols


3. Start now SX-8G V52 software through SX-8G V52 exe.file
4. Press button for Easy Protocol Start screen and load the protocol of interest

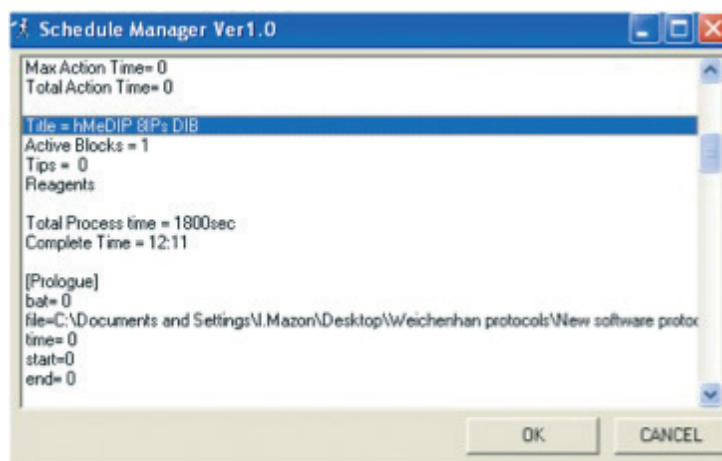

Alternatively, temperature and incubation time for the IP reaction can be modified in an existing protocol by selecting the modify button. The modified protocol can be also saved as new protocol.

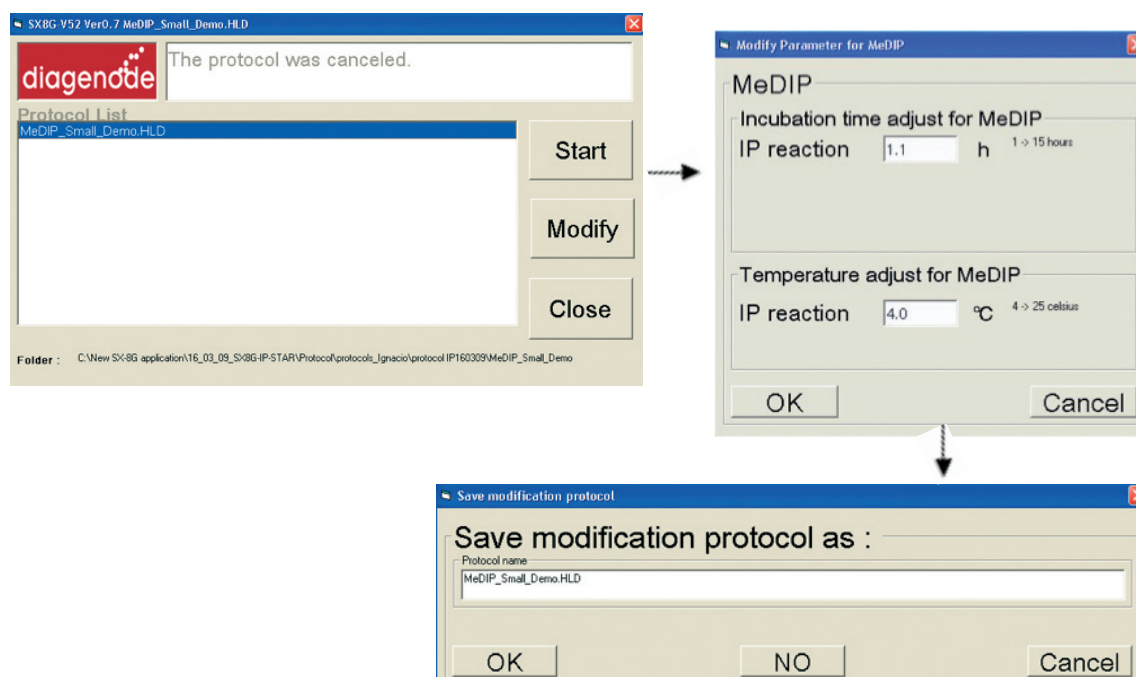

7. The program will run through the following steps: magnetic bead washes, IP and IP washes.

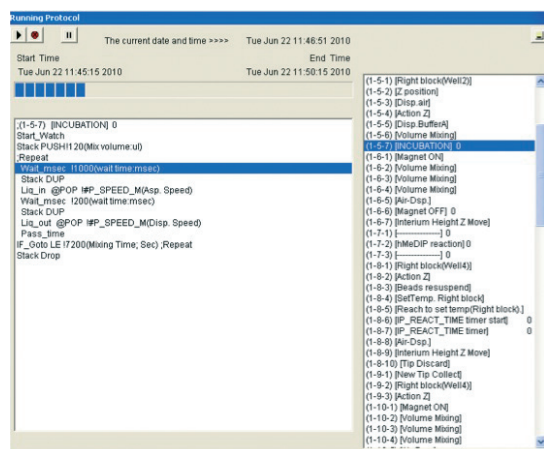


During protocol the next window will be displayed indicating the current protocol step.

8. After the IP washes the following window will be appear.



Follow the next instructions:

1. Add 10 μ l of Input to well 1
2. Add 1 μ l proteinase K to wells 1 and 10
3. Close the tube strip with the corresponding caps
4. Press OK

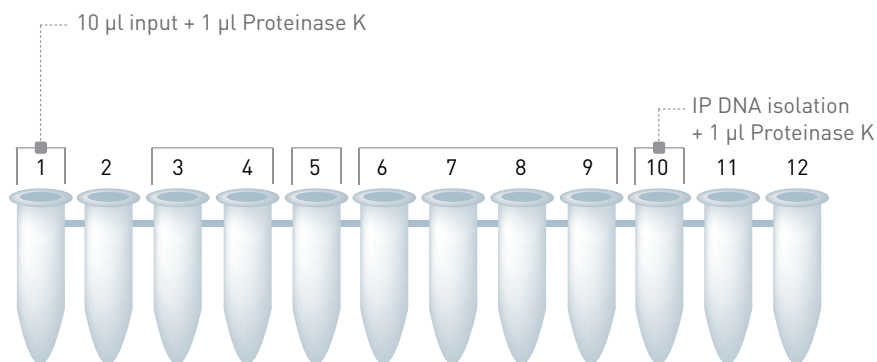

IMPORTANT: Please note that the enriched hydroxymethylated DNA in DIB buffer is single strand DNA that can be directly analyzed by qPCR. For downstream applications such as sequencing or arrays, the enriched methylated DNA needs to be purified by phenol/chloroform extraction and converted to double stranded DNA.

9. The following window will appear:



Close the workstation door and press OK.

The program will move forward to the next steps of the hMeDIP protocol.

10. The SX-8G IP-Star software indicates the end of the protocol.
Collect your immunoprecipitated and isolated DNA
11. Discard magnetic beads by using the DiaMag02 (cat# kch-816-001) or by centrifugation.
12. This is your DNA ready for qPCR.

Shutting down the SX-8G IP-Star

1. Click on File and press End to close the software correctly.
 2. Switch off the computer and its monitor.
 3. Switch off the SX-8G IP-Star Automated System (power switch on the right side).
- Note:** Ensure that the door is closed!

Quantitative PCR & Data Analysis

This last step consists in amplifying and analysing the IP'd DNA.

1. Prepare your qPCR mix using SYBR PCR Green master mix and start out qPCR.

qPCR mix (total volume of 25 µl/reaction:

- 6.50 µl of water
- 12.50 µl of master mix
(e.g.: iQ SYBR Green supermix)
- 1.00 µl of provided primer pair
(stock: 5 µM each: reverse and forward)
- 5.00 µl of isolated DNA or INPUT

	Temperature	Time	Cycle
PCR Amplification	95°C	7 minutes	x1
	95°C	15 seconds	x40
	60°C	60 seconds	
	95°C	1 minute	x1
Melting Curve	65°C and increment of 0.5°C per cycle	1 minute	x60

2. When the PCR is done, analyse the results. Some major advices are given below.

• Data interpretation

The efficiency of hydroxymethyl DNA immunoprecipitation of particular genomic locus can be calculated from qPCR data and reported as a percentage of starting material: % (hmeDNA-IP/ Total input).

$$\% \text{ (hmeDNA-IP/ Total input)} = 2^{[(\text{Ct}(10\% \text{input}) - 3.32) - \text{Ct}(\text{hmeDNA-IP})]} \times 100\%$$

Here 2 is the AE (amplification efficiency); Ct (hmeDNA-IP) and Ct (10%input) are threshold values obtained from exponential phase of qPCR for the hydroxymethyl DNA sample and input sample respectively; the compensatory factor (3.32) is used to take into account the dilution 1:10 of the input. The recovery is the % (hmeDNA-IP/ Total input).

• Background determination


The final goal of IP is to calculate the enrichment in the same IP sample of: 1/ the specific DNA fragments (corresponding to the hydroxymethylated DNA) in comparison with 2/ non-methylated DNA (i.e. negative unDNA control).

• Relative occupancy can be calculated as a ratio of specific signal over background.

$$\text{Occupancy} = \% \text{ input (specific loci)} / \% \text{ input (background loci)}$$

Relative occupancy is then used as a measure of the hydroxymethylation of a specific locus; it provides clues about specificity of the IP. (background loci) corresponds to the signal obtained with one of the unmethylated DNA kit control.

Results

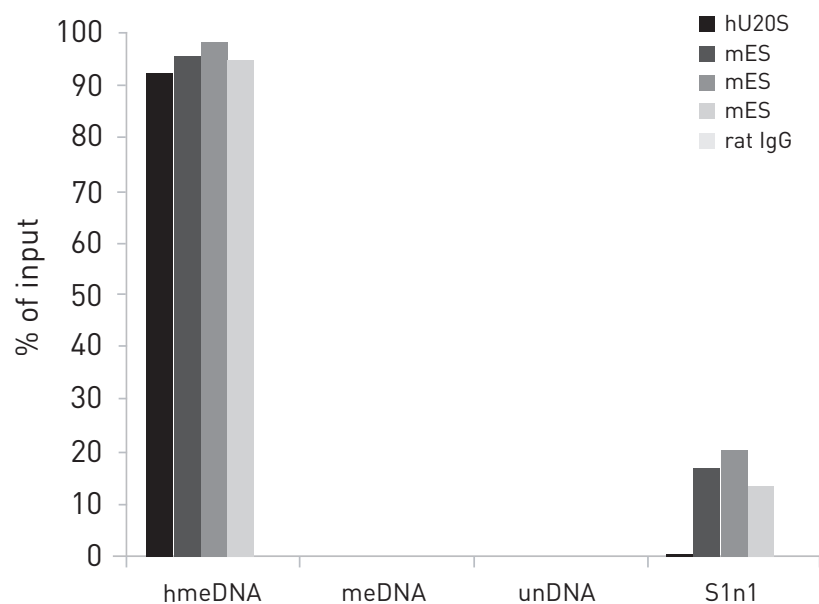

Hydroxymethylated DNA Immunoprecipitation (hMeDIP) was performed using the hMeDIP kit (rat monoclonal antibody) (Cat. No. AF-104-0016) and carried out on the SX-8G IP-Star® Automated System. The Rat IgG isotype antibody was used as negative control. 1 µg of DNA from mouse ES cells was prepared and sonicated with the Bioruptor® to obtain DNA fragments of 300-500 bp. Unmethylated, methylated and hydroxymethylated spike DNA controls (included in the Auto hMeDIP kit) were used in the IP reaction together with the genomic sheared DNA samples. Enrichments were assessed by qPCR using specific primer pairs for the unmethylated, methylated and hydroxymethylated DNA sequences. Sfi1 is a gene that has been identified as being hydroxymethylated using hMeDIP-seq. Human DNA from U2OS was used as negative control DNA.

These results show clearly that hydroxymethylated DNA is specifically immunoprecipitated (hydroxymethylated control and Sfi1) and validate therefore the hMeDIP assay based on Diagenode's hMeDIP kit carried out on the SX-8G IP-Star® Automated System.

Troubleshooting Guide

Error Cause	Remedy
SX-8G IP-Star cannot be switched on	SX-8G IP-Star is not receiving power. Check that the power cord is connected to the workstation and to the wall power outlet.
Computer cannot be switched on	Computer is not receiving power. Check that the power cord is connected to the computer and to the wall power outlet.
SX-8G IP-Star shows no movement when a protocol is started	SX-8G IP-Star is not switched on. Check that the SX-8G IP-Star is switched on.
SX-8G IP-Star shows abnormal movement when a protocol is started	The pipettor head may have lost its home position. In the Software, select "Manual Operation/Home". After confirming that the pipettor head moves to the home position, run the protocol again.
Aspirated liquid drips from the disposable tips	Dripping is acceptable when ethanol is being handled. For other liquids: air is leaking from the syringe pumps. Grease or replace the O-rings. If the problem persists, contact DIAGENODE Technical Services.

Ordering information

Description	Cat. No. (NEW)	Cat. No. (OLD)	Format
SX-8G IP-Star® Compact	B03000002	UH-002-0001	1 unit
Auto True MicroChIP kit	C01010140	/	16 rxns
Auto True MicroChIP & MicroPlex Library Prep Package	C01010141	/	16 ChIP rxns & 12 library prep rxns
MicroPlex Library Preparation kit x12	C05010010	AB-004-0012	12 rxns
Auto Histone ChIP-seq kit protein A x16	C01010020	AB-Auto02-A016	16 rxns
Auto Histone ChIP-seq kit protein A x100	C01010022	AB-Auto02-A100	100 rxns
Auto Histone ChIP-seq kit protein G x16	C01010021	AB-Auto02-G016	16 rxns
Auto Histone ChIP-seq kit protein G x100	C01010023	AB-Auto02-G100	100 rxns
Auto Transcription ChIP kit protein A x16	C01010030	AB-Auto03-A016	16 rxns
Auto Transcription ChIP kit protein A x100	C01010032	AB-Auto03-A100	100 rxns
Auto Transcription ChIP kit protein G x16	C01010031	AB-Auto03-G016	16 rxns
Auto Transcription ChIP kit protein G x100	C01010033	AB-Auto03-G100	100 rxns
Auto ChIP kit protein A x100	C01010011	AB-Auto01-A100	100 rxns
Auto ChIP kit protein G x100	C01010013	AB-Auto01-G100	100 rxns
Auto MeDIP kit x16	C02010011	AF-Auto01-0016	16 rxns
Auto MeDIP kit x100	C02010012	AF-Auto01-0100	100 rxns
Auto hMeDIP kit x16	C02010033	AF-Auto02-0016	16 rxns
Auto MethylCap x48	C02020011	AF-Auto01-0048	48 rxns
Auto IPure kit	C03010010	AL-Auto01-0100	100 rxns

Visit us at one of Diagenode's demo sites or discover our Automated Systems by performing some assays with the help of our R&D and Technical Department.

www.diagenode.com

Diagenode s.a. BELGIUM | EUROPE

LIEGE SCIENCE PARK
Rue Bois Saint-Jean, 3
4102 Seraing - Belgium
Tel: +32 4 364 20 50 | Fax: +32 4 364 20 51
orders@diagenode.com
info@diagenode.com

Diagenode Inc. USA | NORTH AMERICA

400 Morris Avenue, Suite #101
Denville, NJ 07834 - USA
Tel: +1 862 209-4680 | Fax: +1 862 209-4681
orders.na@diagenode.com
info.na@diagenode.com

For a complete listing of Diagenode's international distributors visit:

www.diagenode.com/en/company/distributors.php
For rest of the world, please contact Diagenode s.a.